


## OVERVIEW

### AMS-4085

The Steris 4085 general surgical table handles virtually all the surgical specialties you expect a general surgical table to handle and provides outstanding access for C-arm, permitting high quality images crucial to cardiothoracic, orthopedic and neurosurgical procedures.

- Four-section table top contours to all patient postures
- Powered slide improves C-arm access
- Leg section can be quickly removed for unrestricted perineal access
- Powered, radiolucent kidney elevator enhances lateral procedures
- High weight capacity allows for use with virtually all patients
- Light weight makes the table easy to move
- Impact-resistant, cascading shrouds enable fast, easy cleaning
- Compatible with a wide range of surgical accessories
- Easy-to-use hand control with Return to Level button and Auto Limit Sensor, which prevents unintended section collisions

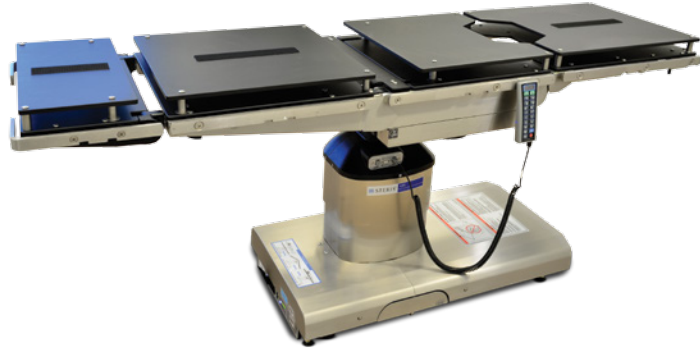

## TECHNICAL SPECIFICATIONS

Length	81" (2057mm)
Width	20" (508mm)
Weight of Operating Table	560 lbs (254 kg)
Maximum Permissible Load (Lift)	1000 lbs (454 kg)
Maximum Permissible Load (Articulation)	600 lbs (272 kg)

## HAND CONTROL ADJUSTMENTS

Height (Including Cushions)	26-45" (660-1143mm)
Lateral Tilt (Left/Right)	20°
Trendelenburg	30°
Reverse Trendelenburg	30°
Top Slide	18" (457.2mm)
Powered Kidney Bridge	4" (150mm)
Back Up/Down	80°/40°
Leg Up/Down	0°/105°
Flex/Reflex	140°/100°
Base Braking (Lock/Unlock)	Yes
Return to Level	Yes

## OPTIONAL ACCESSORIES

Item	Part Number	Description
Accessory Package	ACC-P2	2 Armboards w/2" Pads, 2 Armboard Safety Straps and 1 Standard Patient Restraint Strap
SC Cushion Set	SC-3380	Standard Series - 3 piece Cushion Set
PC Cushion Set	PC-3380	Pressure Management Series - 3 Piece Cushion Set
IGC Cushion Set	IGC-3380	Integra-Gel Series - 3 Piece Cushion Set
Power Cord	PRT-2130	AC Power Cord - 20' Long w/Straight Connector
Hand Pendant	AMS-4085-HC	Remanufactured Hand Control with Cord
X-Ray Tops	RA-5140	X-Ray Panel Set with Stainless Steel Pegs
Perineal Filler Plate	PRT-1230	Surgical Table Perineal Filler Plate



REMANUFACTURED SURGICAL TABLE WARRANTY

2 YEARS PARTS	2 YEARS LABOR	10 YEARS SERVICEABLE GUARANTEE
---------------	---------------	--------------------------------

[www.brytoncorp.com](http://www.brytoncorp.com) • 1-800-567-9500